

Catheter quiver

Standard version and Long version



Optimum protection, hygienic storage

The catheter holder is used for the proper storage of catheters. In order to meet the increasing hygiene requirements on the wards, HERZmed GmbH designed the catheter holder in accordance with the hygiene requirements for sterile supplies (DIN 58953-8: 2010-05).

The stable unit protects against dust, dirt, excessive exposure to light and damage caused by mechanical influences. A dry, sealed and clean interior and smooth, powder-coated surfaces enable optimum hygienic cleaning. The catheter holder and storage box unit is therefore particularly suitable for use in clinics, hospitals, doctors' surgeries and the care sector.

Aluminum, powder-coated in light grey (RAL 7035)
other RAL colors on request



Mounting on standard rail

Wall bracket and accessories optionally available.



Back view



Partition plate



Highlights

- Compact design
- Flexible mounting option
- Transparent lid
- Hygienic storage thanks to closed system
- Flexible height subdivision
- Sealing possible
- Individually adaptable to the customer's needs


Technical details

Transparent lid

The transparent lid provides a quick and easy view into the quiver and allows the filling quantity to be checked. Thanks to the opening angle of 125°, easy filling and removal is guaranteed.



Sealing possible

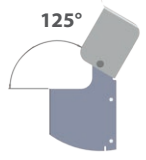

125°

Repainting possible

The catheter holders are supplied as standard in white aluminum RAL 9006 as standard. Repainting in accordance with the RAL color code is possible.


Partition plate inside

A removable divider divides the interior of the cooker diagonally. This division of space makes it possible to separate 2 different catheters.

Conical clips for easy attachment

Various mounting options

The aluminum quiver comes with various mounting options. So that it is always ready for use, it can either be mounted directly on a wall or hung on standard rails.

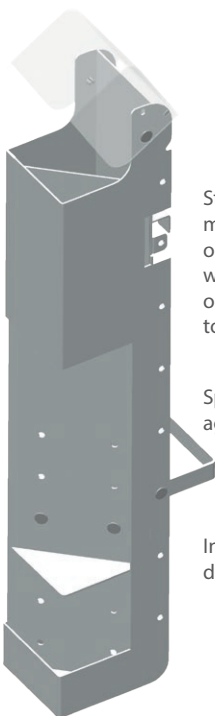

Standard rail bracket for mounting on a standard rail or optionally on the wall with the wall bracket. Other mounting options are possible according to customer requirements.

Spacer bracket for supporting against the wall in two sizes.

Including divider plate with two height dividers for different catheter sizes.

Ergonomic solutions for your workplace

The catheter holder can be individually adapted to customer requirements. Thanks to different interior layouts and various fastening options, it is ready for use anywhere and at any time. It can also be repainted in other RAL colors on request. The product also focuses on ergonomic solutions. This is characterized, among other things, by the optimal removal of the catheters.

Technical specifications

Size and weight Standardversion:

100 x 120 x 703 mm (B x T x H)
2200 g Empty weight

Size and weight Langversion:

100 x 120 x 1052 mm (B x T x H)
3850 g Empty weight

Coloring:

Body: RAL 9006,
Recoating possible on request
Cover: Transparent

Material:

Cover: PETG (glycol-modified polyethylene terephthalate);
Body and parts: White aluminum, powder-coated;
Conical clip: PA66 plastic

Order information

Art.-Nr.	Description
217020	Catheter quiver Standard version, mounted
218000	Catheter quiver Long version, mounted



Body RAL 9006 // **Divider** // **Height divider** Set of 2 incl. 2x Clip // **Lid** incl. lid rivet // **Base** incl. 2x pan-head screws M5 x 10 and 2x nut M5, self-locking // **Standard rail bracket** incl. 2x pan-head screws M5 x 10 and 2x nut M5, self-locking // **Distance holder long and short** (each incl. 2x Clip) // **Instructions** for use and assembly.

Single parts and supplies:

Abb.	Unit	Art. No.	Description	Abb.	Unit	Art. No.	Description
	Piece	227027	Catheter quiver Distance holder short (incl. 2x clip)		Piece	227023	Catheter quiver height divider set of 2 incl. 2x clip
	Piece	227028	Catheter quiver Wall mount VESA 4x pan-head screws M5 x 10 and 4x nut M5, self-locking)		Piece	227024	Catheter quiver lid
	Piece	227029	Cone clip		Piece	227025	Catheter quiver Standard rail bracket incl. 2x pan-head screws M5 x 10 and 2x nut M5, self-locking
	Piece	029.0446	Universal claw (for wall mount)		Piece	227026	Catheter quiver Distance holder long (incl. 2x clip)
	Piece	227022	Catheter quiver Standard version divider	Art.-No. Description			
				029.0446			Universal claw (for wall mount)
				227031			Conical clip (spare parts kit) consisting of 10 pieces
				227032			Screw set made of VA (spare parts kit), 10x pan-head screws M5 x 10 and 10x nut M5, self-locking

Vertrieb/Sales Office:

HERZmed GmbH
Rottenbucher Straße 9
82497 Unterammergau

Tel.: +49 (0) 8822 824 0002
Fax: +49 (0) 8822 824 0006

Verwaltung/Managment:

HERZmed GmbH
Herkulan-Schwaiger-Gasse 16
82487 Oberammergau

Mail: herzmed@herzmed.de
Web: www.herzmed.com

The products shown are only an excerpt of our product range. Exact details and product information are available on request or at www.herzmed.de. All prices marked with an asterisk are exclusive of VAT. All brand names used remain the property of the respective brand owner. We assume no liability for errors, possible misprints and color deviations. Our general terms and conditions apply, which you can also view at www.herzmed.de.

